Jordan CLEAN

Electrical Toothbrush

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Valmistaja: **WILFA AS, Industriveien 25, NO-1481 Hagan, Norja**

JORDAN CLEAN ELECTRIC TOOTHBRUSH

The Jordan Clean series is a new generation of electric toothbrushes that through SoftClean™ Technology delivers the right combination of oral hygiene features. Inspired by a Scandinavian design approach for simplicity and user-friendliness, the series gives an effective clean from an electric toothbrush and the ease of use of a manual one.

SAFETY

Please note the following instructions when using the appliance:

- This appliance is intended for domestic use only.
- This appliance is designed solely for cleaning and caring for human teeth. All other uses are expressly prohibited.
- Do not use if the appliance or attachments is visibly damaged.
- When the appliance is being charged, on no account should this or the charging unit come into contact with water or any other liquid.
- Keep the appliance out of direct sunlight, excessive heat and open flames.
- Do not drop the appliance.
- Keep the appliance out of the reach of children.
- Do not open the appliance under any circumstances.
- No warranty claims are accepted for damage caused by improper handling.
- Cleaning and user maintenance shall not be made by children without supervision.
- This appliance can be used by children aged from 8 years and above. And persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- The Electric Toothbrush can only be charged by the charging unit provided with the appliance.

SPECIAL FEATURES:

- **Softclean™ technology** gives you a unique combination of highly effective bristles, and a smaller head for effective cleaning in an effective way. Each head is divided into three cleaning zones with specialized efficiencies per zone.
- **High oscillating speed** at about 8800 rpm for effective cleaning performance.
- **2-min timer** with click signal for each 30 sec brushing time for reminding you to move to different teeth area. It will click 2 times at 30/60/90 sec. and click 4 times at 2 minutes (recommended brushing time). The brush will stop automatically after 10 minutes for protecting your teeth and preserve battery capacity by accidental switch-on when travelling.
- **Low power indication** when battery capacity goes down to about 20%, the indicator blinks quickly to remind you to recharge the brush in

order to keep full power.

- **Built-in high-capacity nickel-metal hydride battery** for long use. The toothbrush can operate for 60 minutes, which is equivalent to 2 weeks if you brush twice a day and spend 2 minutes brushing each time.
- **Over-charge and over-discharge Protection** circuit is equipped to preserve the battery life time. When fully charged the toothbrush will automatically cut off. When running low on battery the toothbrush will automatically stop.

SPECIFICATIONS

Charging unit: 100-240V~ 50/60Hz 1W.
Rechargeable power unit is 2.4Volts. Maximum operation time is approximately 60 minutes/2 weeks when fully charged. When rechargeable batteries are running low, the indicator will blink to tell you to recharge the unit immediately. If you continue to use the unit for some time without recharging it, the unit may stop automatically. This is our battery protection feature to preserve the batteries life time. You need to recharge immediately.

CHARGING THE UNIT

The charging unit can be connected to any voltage outlet - even to shaver sockets - between 100-240 volts AC (~). The indicator will blink with light/shade (breathing mode) alternately, when the power unit is inserted to charging base which is plugged-in. It will light up WHITE when batteries are fully charged (18 hours). The power unit may warm up somewhat during recharging. This has no significance on the unit's operation. Over-charging is impossible.

PRESERVING THE RECHARGEABLE BATTERIES


To maximize the capacity and life of the rechargeable batteries, they should be fully charged before first use and then fully discharged through use. Then recharge to full capacity again. This discharging/charging procedure should be done every six months. The Electric Toothbrush also has built-in LOW BATTERY INDICATION feature. If the battery power often runs too low this will shorten the battery life. The indicator will blink to remind you to charge. If you don't charge another "Automatic Shut off system" will stop the brush. You need to fully recharge before using again.

BEFORE USE

Before using your Jordan Clean toothbrush for the first time, charge the unit full. After the power unit has been fully charged, attach the brush head to the power unit. The brush head must face the same direction as the switch after twist in. Jordan Clean has only one switch by pressing for the desired functions. In case you put the power unit into the charging base when it is working, it will automatically stop.

REMOVAL OF RECHARGEABLE BATTERY

This toothbrush is using Ni-Mh rechargeable battery. When the toothbrush reaches the end of its

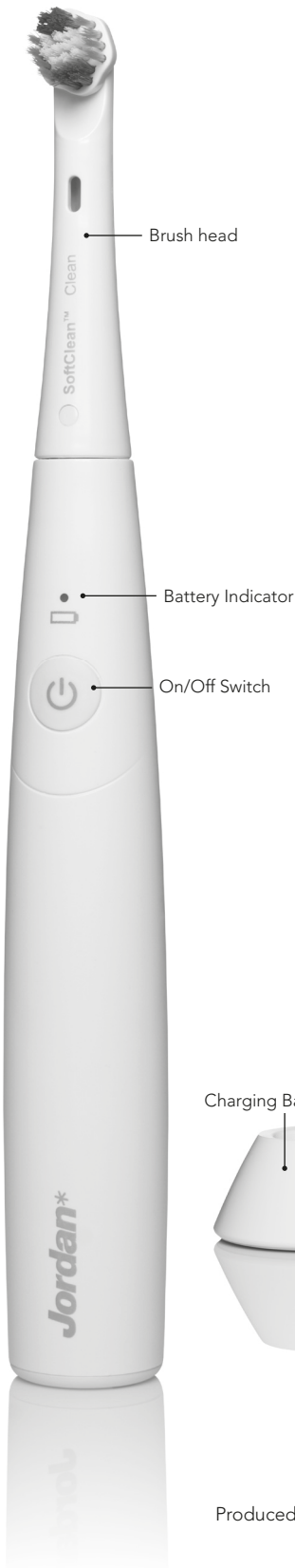service life, the battery must be removed from the appliance before it is discarded. The appliance must be disconnected from the mainpower supply before removing the battery. The battery must be disposed of safely. Use a flathead screwdriver to remove the internal body, then remove the battery and dispose of it via a return and collection system or contact the retailer from whom the product was purchased. They operate a take-back system for safe recycling.

PROTECTING THE ENVIRONMENT

Correct disposal of this product.

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

POWER SUPPLY

Adaptor:

100-240V~, 50/60Hz, 1W.

Appliance:

Nickel metal hydride battery 2.4V.

Produced by: **WILFA AS, Industriveien 25, 1481 Hagan, Norway.**

Markedsføres i Norge av:
WILFA AS
Industriveien 25
Postboks146, Hagan
1481 Hagan, Norway
Tlf: 67 06 33 00
support@wilfa.com
www.wilfa.com

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WILFA AB
Oy Wilfa Suomi Ab